


## USER'S GUIDE

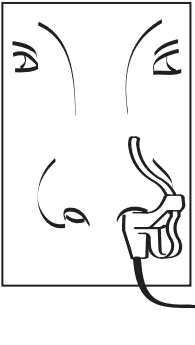

# Vielight® The Life Light® 810 Infrared

## OPERATING INSTRUCTIONS

### Treatment Sessions

Each treatment session lasts 25 minutes. One session per day is recommended. It can be used at anytime but bedtime may be the most convenient. The unit will automatically switch off at the end of each session.

### Using Vielight 810 Infrared

1. Open the back cover by pressing the release button (h), hold it down, and simultaneously slide open the back cover (j). Insert a 1.5V AA battery.
2. Fully insert the jack plug (f) into the socket (g).
3. Set the applicator lens (d) deep inside nose, with the clip (e) on the outside.
4. Press the button (b) to switch the unit on. A single beep will sound.
5. The unit will remain switched on for 25 minutes. Indicator lights (c) will blink throughout the session. The number of blinking indicator lights will reduce over time, which is indicative of the remaining time.
6. At the end of the session, the unit will beep three times and then switches off. The unit may be switched off at mid-session by pressing button (b).

### Replacing the Battery



1. Replace the battery when the Control unit beeps 5 times continuously, or when the low battery indicator blinks.

Note: Remove the battery if the equipment is not likely to be used for some time.

2. To access the battery compartment, open the back cover by pressing the release button (h), hold it down, and simultaneously slide open the back cover (j). Insert a 1.5V AA battery.

### Care & Cleaning

1. To clean the Applicator, gently wipe the surface with a cotton swab or a soft cloth moistened with hydrogen peroxide. You can wipe the applicator regularly with clean water.
2. Use a soft, dry cloth or cotton swab to clean the other surfaces. Do not use abrasive cleaners. Do not submerge the Applicator or Control Unit in water.
3. Store the equipment in a dry location, free from dust, away from direct sunlight. The storage temperature should be within 10°C to 40°C.
4. Clips (e) and (j) are made of tough, high-impact plastic materials. However, excessive bending can cause permanent deformation or breakage.

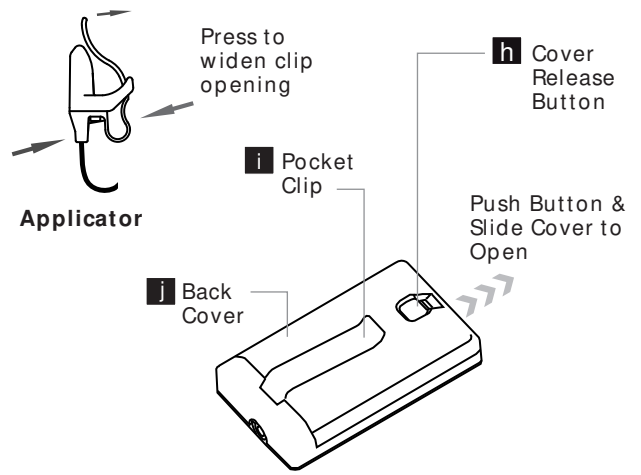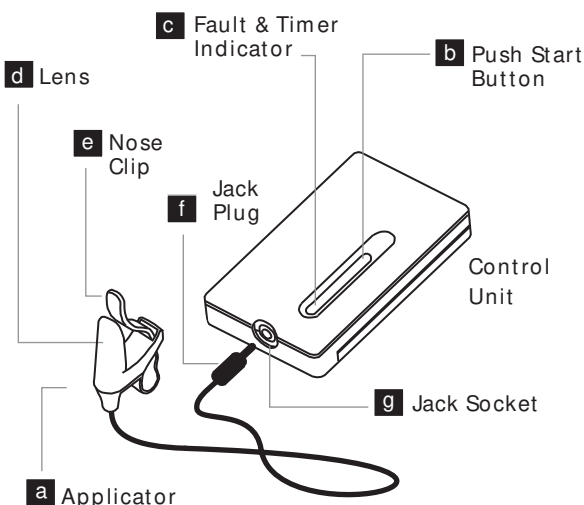

### Indicator Legend

	Battery low	Replace the battery
	The light source is not functioning	Check that the applicator is completely plugged in. If indicator persists, replace the applicator.

### Output Visual Check

Whilst the light beam is invisible under normal light, it is visible as a dim flickering red light in darkness.

### Removing Battery Cover

1. Hold the casing as illustrated, and press the cover release button.
2. While holding down the release button, push to slide open the cover.
3. Insert the battery as shown

### Product Classification

- Internally powered equipment
- Type B applied part - nose applicator
- Ordinary equipment: IPX0
- Equipment is not suitable for use in the presence of a flammable anesthetic mixture with air or with oxygen or nitrous oxide
- No significant electromagnetic or other interference between Vielight 810 Infrared and other equipment

### Disposal

To protect the environment, dispose of used batteries at appropriate collection sites according to national or local regulations.

Plastic parts are resistant to biodegradation and should be recycled appropriately.



## Technical Specification

Light Source:	Solid-state, 810 nm light-emitting diode
Power:	Regulated, pulsed at 10 Hz
Auto Off Timer:	25 minutes
Indicators:	Battery low Applicator fault Application timer indicator
Battery Type:	1 x 1.5V AA
Battery Life:	Approx. 10 hours with alkaline battery

Product by Vielight Inc., Canada



Follow operating instructions  
Symbol for advising operator that operating instruction must be followed



Type B applied part  
Symbol for protection against electric shock for nose applicator

## Limited One-Year Warranty

Your Vielight 810 Infrared is warranted to be free of defect in material and workmanship under normal use.

This warranty is limited for a period of one year from its original purchase date in which repairs or exchange of parts according to what is deemed necessary by the company will be made free of charge.

This warranty does not cover damage by the following;

- Misuse, abuse or alteration to the product
- Damage caused by accident or neglect
- Servicing or repairs made by unauthorized persons

WE SHALL NOT BE RESPONSIBLE FOR ANY INCIDENTAL, SPECIAL, OR CONSEQUENTIAL DAMAGES RESULTING FROM THE USE OF THIS PRODUCT.

V810IR310725